




Alpha-208

Double Wheels Grinder & Polisher

GRINDING/POLISHING

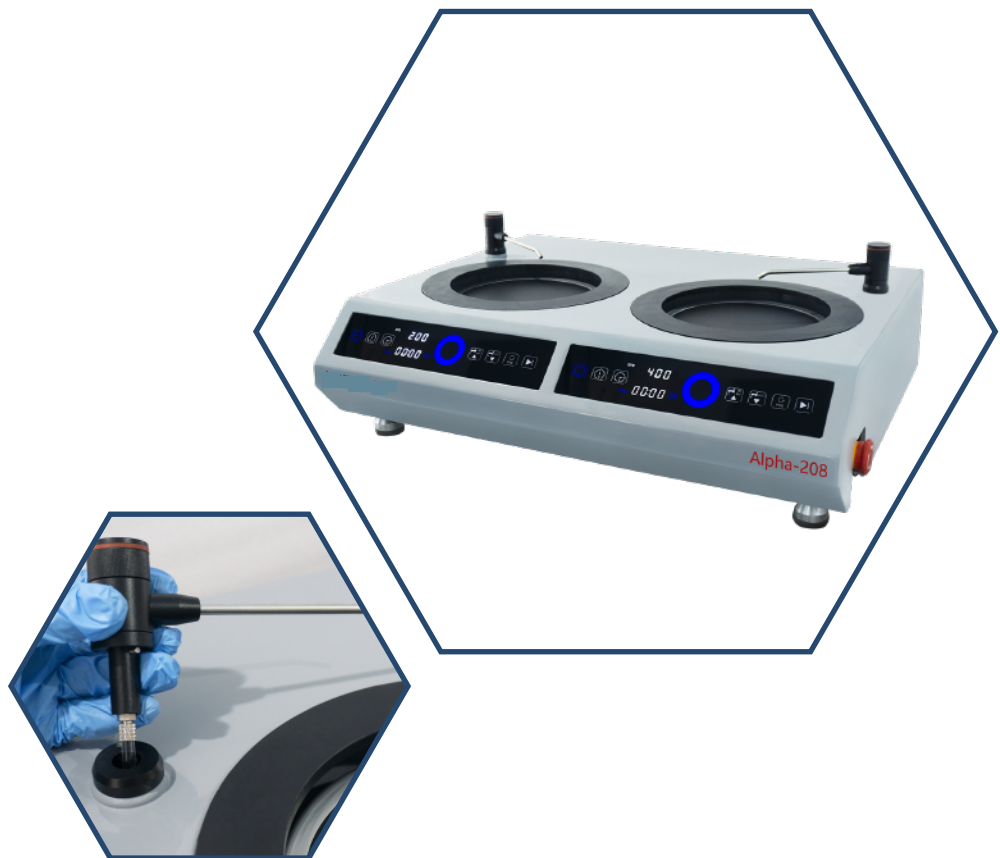

2022

04 / 15

Alpha-208 Double Wheels Grinder & Polisher

Alpha-208 is a newly designed grinder & polisher, that have individual control panel and double wheels of 8"(203mm) diameter. This machine is compact and easy to operate. Adopted the current popular touch key-press panel, and gesture-sensing function, stable in performance, especially suitable for grinding and polishing all kinds of samples.




Features:

- ◆ 750W servo motor to maintain high and constant torque output even at low speed;
- ◆ Glass touch press panel, easy to use and simple to operate;
- ◆ Gesture-sensor function makes operate easily;
- ◆ The spindle adopts fully sealed design, even if the water has no time to discharge from the drain hose, it will not leak into the equipment;
- ◆ Retractable water spout for convenient cleaning and maintain;
- ◆ The disc is easy to remove without any tools;
- ◆ ABS housing, more strength and no rust.



Control Panel Introduction:

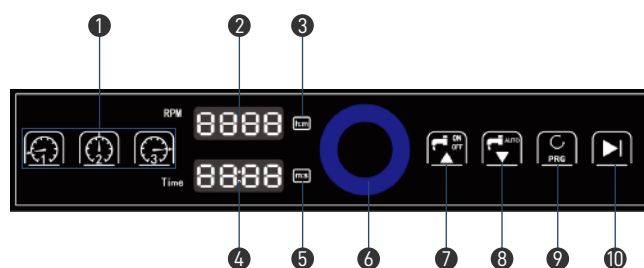

- | | |
|--------------------------------|----------------|
| ① Preset Speeds | ② Display RPM |
| ③ Timer format:Hours:Minutes | ④ Timer |
| ⑤ Timer format:Minutes:Seconds | ⑥ Speed Adjust |
| ⑦ Water ON/OFF | |
| ⑧ Set water to Auto mode | ⑨ CW/CCW |
| ⑩ Start/Stop | |


Technical Specifications:

Model		Alpha-208
Wheel	Diameter	8"(203mm)
	Speed	10-1000rpm
	Direction	CW/CCW
	Torque	10N.m
Electronic	Voltage/Frequency	220V/50Hz (1Ph/N/PE)
	Rated Power	1500W
	Rated Current	8A
Control	Panel	Glass & Touch
	Preset Speeds	Yes
	Timer	Yes
	Gesture -Sensor	Yes
Supply	Max. water inlet	<0.3MPa
	Inlet type	G 3/4 (Male)
	Drain Hose	OD=Φ32mm(1-1/4 inch)
	Recirculating Tank	Optional
Noisy	Measure at 1 meter away	<50dBA
Operating Environment	Temperature	-20°C - 60°C (-4°F -140°F)
	Humidity	0-95%RH
Dimension	WxDxH	672×575×280mm
Weight		34kg




Dimension Figure

Unite: mm

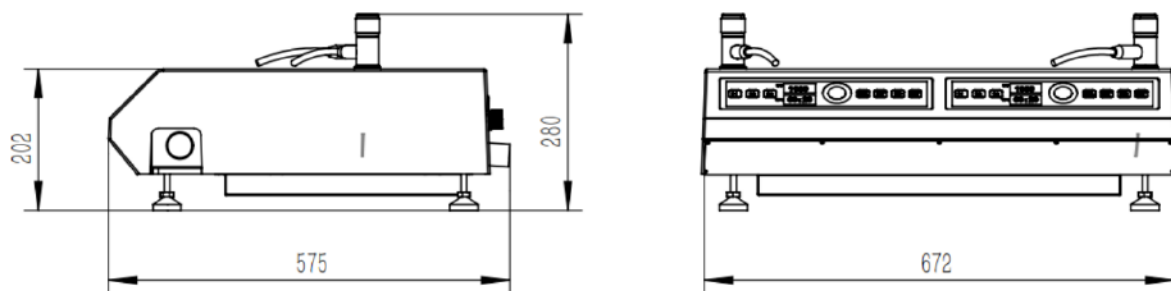


Related Consumables & Accessories:

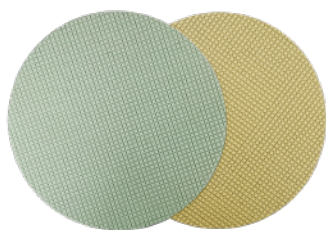

DiaRe Diamond Disc



Abrasive Paper

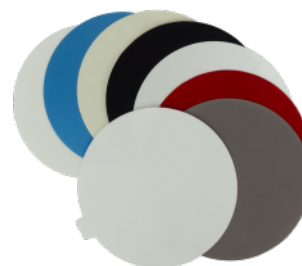

Polishing Cloth



Diamond Suspension



Alumina/Silica suspension



Adapter Plate System



Purchase list:

DiaNi Diamond Disc:

	8in DiaNi (200mm)	10in DiaNi (250mm)	12in DiaNi (300mm)	Package
80	04.01.110S	04.01.210S	04.01.310S	1PC
120	04.01.120S	04.01.220S	04.01.320S	1PC
200	04.01.125S	04.01.225S	04.01.325S	1PC
400	04.01.130S	04.01.230S	04.01.330S	1PC
800	04.01.140S	04.01.240S	04.01.340S	1PC
1200	04.01.150S	04.01.250S	04.01.350S	1PC

DiaRe Diamond Disc:

	8in DiaRe (200mm)	10in DiaRe (250mm)	12in DiaRe (300mm)	Package
80	04.02.110S	04.02.210S	04.02.310S	1PC
120	04.02.120S	04.02.220S	04.02.320S	1PC
200	04.02.125S	04.02.225S	04.02.325S	1PC
400	04.02.130S	04.02.230S	04.02.330S	1PC
800	04.02.140S	04.02.240S	04.02.340S	1PC
1500	04.02.150S	04.02.250S	04.02.350S	1PC
3000	04.02.160S	04.02.260S	04.02.360S	1PC
5000	04.02.170S	04.02.270S	04.02.370S	1PC

Silicon Carbide Disc:

	8in PSA (200mm)	10in PSA (250mm)	12in PSA (300mm)	8in Plain Back (200mm)	10in Plain Back (250mm)	12in Plain Back (300mm)	Package
120	03.04.203	03.04.503	03.04.803	03.03.103	03.03.403	03.03.703	50PCS
180	03.04.204	03.04.504	03.04.804	03.03.104	03.03.404	03.03.704	100PCS
400	03.04.207	03.04.507	03.04.807	03.03.107	03.03.407	03.03.707	100PCS
800	03.04.209	03.04.509	03.04.809	03.03.109	03.03.409	03.03.709	100PCS
1200	03.04.211	03.04.511	03.04.811	03.03.111	03.03.411	03.03.711	100PCS
2000	03.04.213	03.04.513	03.04.813	03.03.113	03.03.413	03.03.713	100PCS
2500	03.04.214	03.04.514	03.04.814	03.03.114	03.03.414	03.03.714	100PCS
2400	03.04.215	03.04.515	03.04.815	03.03.115	03.03.415	03.03.715	100PCS
4000	03.04.216	03.04.516	03.04.816	03.03.116	03.03.416	03.03.716	100PCS

Alumina Disc:

	8in PSA (200mm)	10in PSA (250mm)	12in PSA (300mm)	8in Plain Back (200mm)	8in Plain Back (250mm)	8in Plain Back (300mm)	Package
120	03.02.403	03.02.503	03.02.803	03.01.103	03.01.103	03.01.703	50PCS
180	03.02.404	03.02.504	03.02.804	03.01.104	03.01.104	03.01.704	100PCS
320	03.02.405	03.02.505	03.02.805	03.01.105	03.01.105	03.01.705	100PCS
600	03.02.407	03.02.507	03.02.807	03.01.107	03.01.107	03.01.707	100PCS
800	03.02.408	03.02.508	03.02.808	03.01.108	03.01.108	03.01.708	100PCS
1200	03.02.410	03.02.510	03.02.810	03.01.110	03.01.110	03.01.710	100PCS

Adapter Plate System:

	AD-PSA Plate	AD-Foil Plate	AD-PB Plate	Magnetic Plate	Package
200mm	04.08.410	04.08.422	04.08.310	04.08.010	1PC
250mm	04.08.413	04.08.430	04.08.320	04.08.020	1PC
300mm	04.08.420	04.08.432	04.08.330	04.08.030	1PC



Purchase list:

Polishing Cloth:

		8in (203mm)	10in (254mm)	12in (305mm)	Package
Rough Polishing	YS-JP White Fiber Polishing Cloth(PSA)	05.07.710	05.07.720	05.07.730	10pcs
	GF-JP White Synthetin Silk Polishing Cloth (PSA)	05.10.010	05.10.020	05.10.030	10pcs
Fine Polishing	SR-ZP Red Synthetic Fabric Polishing Cloth (PSA)	05.02.150	05.02.160	05.02.170	10pcs
	SC-JP White Natural Silk Polishing Cloth (PSA)	05.05.510	05.05.520	05.05.530	10pcs
	CS-JP Blue Synthetin Silk Polishing Cloth (PSA)	05.12.110	05.12.120	05.12.130	10pcs
Final Polishing	ZN-ZP Black Polyurethane Polishing Cloth (PSA)	05.01.110	05.01.120	05.01.130	10pcs
	ET-JP White Synthetic Velvet Polishing Cloth (PSA)	05.03.210	05.03.220	05.03.230	10pcs
	YR-JP Brown Synthetic Rayon Polishing Cloth (PSA)	05.04.310	05.04.320	05.04.330	10pcs
	Final Red Final Precision Polishing Cloth(Adhesive Back)	05.13.110	05.13.120	05.13.130	10pcs

Diamond polishing suspension and paste:

Diamond Suspension		Micrometer	Diamond Paste		Diamond Spray
Monocrystalline (500ml)	Polycrystalline (500ml)		Monocrystalline (5g)	Polycrystalline (10g)	Monocrystalline (350ml)
06.01.110	06.02.210	0.25um	–	–	–
06.01.120	06.02.220	0.5um	–	–	06.03.115
06.01.130	06.02.230	1um	06.10.030(Blue)	06.10.110(Black)	06.03.320
–	–	2.5um	–	–	06.03.325
06.01.140	06.02.240	3um	06.10.040(Green)	06.10.130(Black)	–
–	–	3.5um	–	–	06.03.330
–	–	5um	–	–	06.03.335
06.01.150	06.02.250	6um	06.10.050(Yellow)	06.10.160(Black)	–
06.01.160	06.02.260	9um	06.10.060(Red)	06.10.190(Black)	–
–	–	10um	–	–	06.03.340

Precision polishing suspension and polishing powder:

	0.05µm	0.3µm	1µm
Alumina Suspension(500ml)	06.04.210	06.04.220	06.04.230
Alumina Polishing Powder(1LB)	06.13.011	06.13.021	06.13.031
	500ml	1000ml	5000ml
SO T401 Non-Crystallizing Silica(50nm)	06.05.310	06.05.320	06.05.330
100C30 silica suspension(50nm)	06.08.310	06.08.311	06.08.312
SO A539 silica suspension(50nm)	06.07.520	06.07.523	

Polishing lubricants:

PL-W Polishing Lubricants	06.09.510	06.09.520	500ml/1000ml
---------------------------	-----------	-----------	--------------



Útiles y Máquinas Industriales, s. a.

Pol. Ind. Ugaldeguren I, Parc. P 3-II, Pab. 7 - 48170 ZAMUDIO - CIF : A48010821

DELEGACIÓN BIZKAIA - 94 446 62 50 - info@umi.es

DELEGACIÓN MADRID - 91 678 46 48 - madrid@umi.es